



CSBio

CS136X Research Scale Peptide Synthesizer



INNOVATIVE.

RELIABLE.

TRUSTED.

CSBio is a leading peptide and peptide synthesizer manufacturing company located in Silicon Valley in Menlo Park, California.

For over 25 years, CSBio has been providing automated peptide synthesizers and cGMP peptides to the global pharmaceutical community.

CSBio builds peptide synthesizers from research scale to commercial scale manufacturing and utilizes all these instruments in-house. Our company is the only peptide synthesizer manufacturer that also manufacturers peptides from mg to kg.

The CS136X is the most popular and purchased instrument in our portfolio. Customers routinely praise the individual customization of their system that perfectly fits their needs.





The engineering team at CSBio takes meticulous care to ensure our instruments are designed and manufactured to the highest specifications. We strive to uphold this standard across all designs, as well as in the selection of the highest quality components and materials for our products.

The CS136X is our most robust and flexible research scale peptide synthesizer, used for R&D, process development, and small-scale manufacturing. The instrument can be found in large pharmaceutical companies and universities across the globe. With all CSBio systems built to order, our synthesizers can be customized to the user's specific needs.

Reliability

CSBio has more synthesizers running over 10 years in the field than any other manufacturer; with 8 of the top 10 pharma companies in the world utilizing CSBio instruments. Many of our instruments in the field have been operational for 15+ years requiring little to no service. Our instruments are constructed with high quality components, feature peptide knowledge enabled integration and are designed for ease of use, as well as ease of repair.

Competitively Priced

Simplicity in design is difficult to achieve. With over 25 years of experience in honing the CS136X, CSBio has been able to achieve a simple and economical, yet extremely effective and reliable design. This allows us to pass these savings directly down to our customers, so they can start up their lab at a lower cost and keep it running for years. Unlike competitors, we prioritize our efforts on the products, not on marketing fluff.

Fast Synthesis with Intuitive Software

Our software and user interface is the most flexible on the market and has been optimized by our in-house chemists who use our instruments daily to synthesize peptides in our R&D and cGMP facilities. CSBio peptide synthesizers operate using the same software across all peptide synthesizers scale up to pilot and commercial scale easy.



Designed with Chemists in Mind

Complete Instrument Synthesis Control

Temperature Control up to 90°C | Using jacketed RV's the 136X employs uniform heating for controlled temperatures with no localized spikes or racemization. Circulating heated water around the reaction vessel during the entire synthesis, not just during coupling. This keeps the growing peptide chain extended, allowing for fast and high purity synthesis. Our system integrated circulator can provide consistent temperatures up to 90°C with the ability to turn on/off the circulator for specific AAs in the sequence.

Adjustable Scale | The 136X has a synthesis range of 0.1 – 4mmol. RV sizes can range from 20 – 250mL, translating to about 100mg - 7g of starting resin.

Flexible Chemistry | This system allows for unlimited automated couplings and is suitable for any type of chemistry, including Fmoc or TBoc, special amino acids as well as multiple activators including HBTU/DIEA and HoBT/DIC, or any other type of solid phase chemistry.

Superior Mixing | Our research scale instruments feature patented 180° inversion mixing with nitrogen bubbling. This type of mixing is the best on the market providing more resin to solvent interaction than vortex mixing, oscillation mixing, or microwave. Delivering high quality crude purities without the problems associated with microwave heating.

High-Speed Synthesis


CSBio peptide synthesizers can consistently produce high purity peptides with time-efficient protocols.

Peptide Synthesis Capabilities Every 30 hours

	CS136X		CS136M
Reaction Vessels	1 RV	2RV	6 RV
Total Peptides Synthesized (20-mer, 0.75 mmol scale)	5	10	18

Flexible User-Friendly Automated Synthesis

CSBio peptide synthesizers run on our in-house developed and proprietary CSPEPM software that is naturally intuitive for users, developed by chemists for chemists. Anything that can be done with manual synthesis, can be automated with our peptide synthesizers.

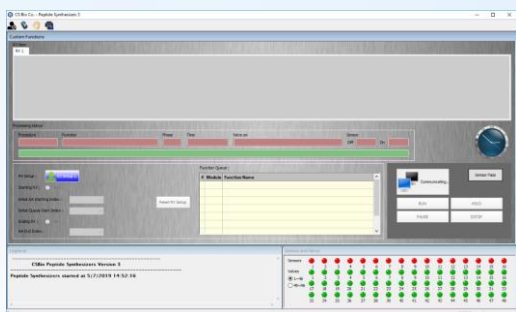

Our software gives the user the ability to:

- Synthesize unlimited automated couplings.
- Build their own step by step protocol or use one of our standard protocols.
- Couple specific AA with different activators in the same sequence.
- The option to cap or double/triple couple individual AA in the sequence.

Custom Options

High Throughput; 3 synthesizers in 1 | The 136X can have up to 3 functional RVs that can make multiple peptides at a time. This can double the output of the instrument. The multiple RVs can simultaneously run peptides that are different sequences, different scales, or even operate under different protocols.

Increased Capacity | We can equip your synthesizer with additional amino acid or solvent bottles. We can also adjust the size of the solvent bottles upon request.

UV Monitoring | Constant detritylation detection to verify deprotection of the amino acids.

GMP Compatible | An upgraded version of the software is compliant to 21 CFR part 11 for GMP manufacturing. As required, we also provide Installation/Operation Qualification (IQ/OQ) for the system.



PTFE Teflon Premium grade virgin PTFE, USDA 21 CFR 177.1550 Certified





GET IN
TOUCH



CSBio Instrumentation Co.
20 Kelly Court
Menlo Park, CA, 94025
U.S.A



+1 650 322 1111
+1 650 322 2278



instrument@csbio.com



www.csbio.com